



671 Extra Ergolift

Powered stair chair



FEATURES

[Visit Extra Ergolift website](#)

Extra Ergolift is the powered version of Extra track assisted stair chair with the same main structural features: lightweight aluminium alloy frame, rear bar and front telescopic handles, practical covered armrests and frontal turnable footrest, long and high-grip sliding tracks which always allow to work on three steps and 4 wheels, two rear fixed Ø 200 mm and two front rotating Ø 125 mm diameter. The real innovation of Extra Ergolift is given by the electric motor that moves the tracks up and down the stairs without any effort by the operator. The motor is activated by a wireless console. The battery ensures an high autonomy and, in case of sudden discharging, this device will allow to safely completing the service. **Extra Ergolift loads up to 230 kg.** It can be used under rain or adverse weather conditions.

Weight 32 kg.

28V battery

Autonomy: 18 ramps by 20 steps with a load of 130 kg

Charging time: max 1 hour



VIDEO MeBer EXTRA ERGOLIFT

VIDEO TUTORIAL

MAIN FEATURES

Cardio	
Max length (cm)	150
Length (cm)	84
Width (cm)	55
Height (cm)	104
Epaisseur pliée (cm)	29
Weight (kg)	32
Loading capacity (kg)	230
Material	Aluminium alloy

IMAGES



SEE ALSO



679
Narrow foldable sedan chair



Extra 670/BR
Track assisted stair chair with arms



Extra Ergolift FX 675
Fixed frame powered Stair Chair
suitable for transportation

ACCESSORIES



612/AN-2/MEB

Adjustable thoracic abdominal belt orange/black with MeBer band and 4 points fastening for self loading stretchers and sedan chairs



612/N-2/MEB

Adjustable thoracic abdominal belt black with MeBer band with 4 points fastening for self loading stretchers and sedan chairs



975

"Ergolift fixing system" - Fixing system for chair Extra Ergolift



12098

Black back rest and seat cushion for Extra chair



671-001

Transpire, dust proof black nylon cover for Extra Ergolift chair




12016-S

Phleboclysis stand (Extra chairs)



670-001

Orange/black head immobilizer for Extra stair chair

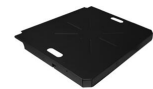

930L

Seat in black embossed painted aluminum for 670/BR and 671 extra chairs



056L

Back in black embossed painted aluminum for 670/BR and 671 extra chairs



670-004

Castor with brake diam. 200 mm for Extra chair



566

Black bag with 3 belts in 2 pieces made of black ribbon with orange MeBer logo for Extra chairs



171J

Removable left lever with ergonomic handle for Extra chairs.

STANDARD ACCESSORIES



12092

Battery MILWAUKEE M28 RED 5,0 Ampere for Extra Ergolift chair



12090

Battery charger MILWAUKEE M28C for rechargeable battery M28 RED Extra Ergolift chair



566

Black bag with 3 belts in 2 pieces made of black ribbon with orange MeBer logo for Extra chairs